


Multilocus PCR-RFLP typing of *Toxoplasma gondii*

Chunlei Su Lab, Microbiology, University of Tennessee, Knoxville

Reference strains


		5'+3'	alt.									
	SAG1	SAG2	SAG2	SAG3	BTUB	GRA6	c22-8	c29-2	L358	PK1	Apico	
GT1,RH88	I	I	I	I	I	I	I	I	I	I	I	
PTG	II or III	II	II	II	II	II	II	II	II	II	II	
CTG	II or III	III	III	III	III	III	III	III	III	III	III	
TgCgCa1 (cougar)	I	II	II	III	II	II	II	u-1	I	u-2	I	
MAS	u-1	I	II	III	III	III	u-1	I	I	III	I	
TgCatBr5	I	III	III	III	III	III	I	I	I	u-1	I	

SAG1


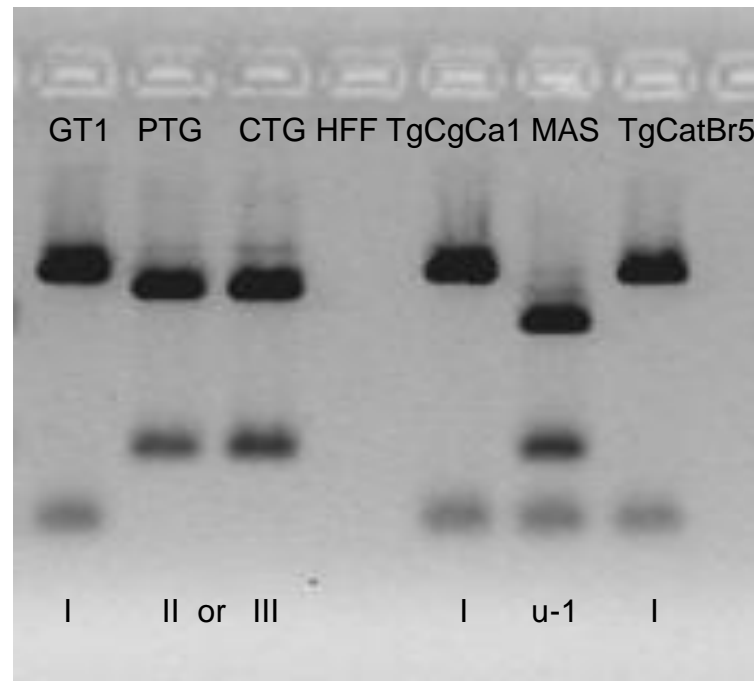


Sau96I + HaeII (double digest) NEB4, BSA, 37C 1 h, 2.5% gel

(5'+3') SAG2

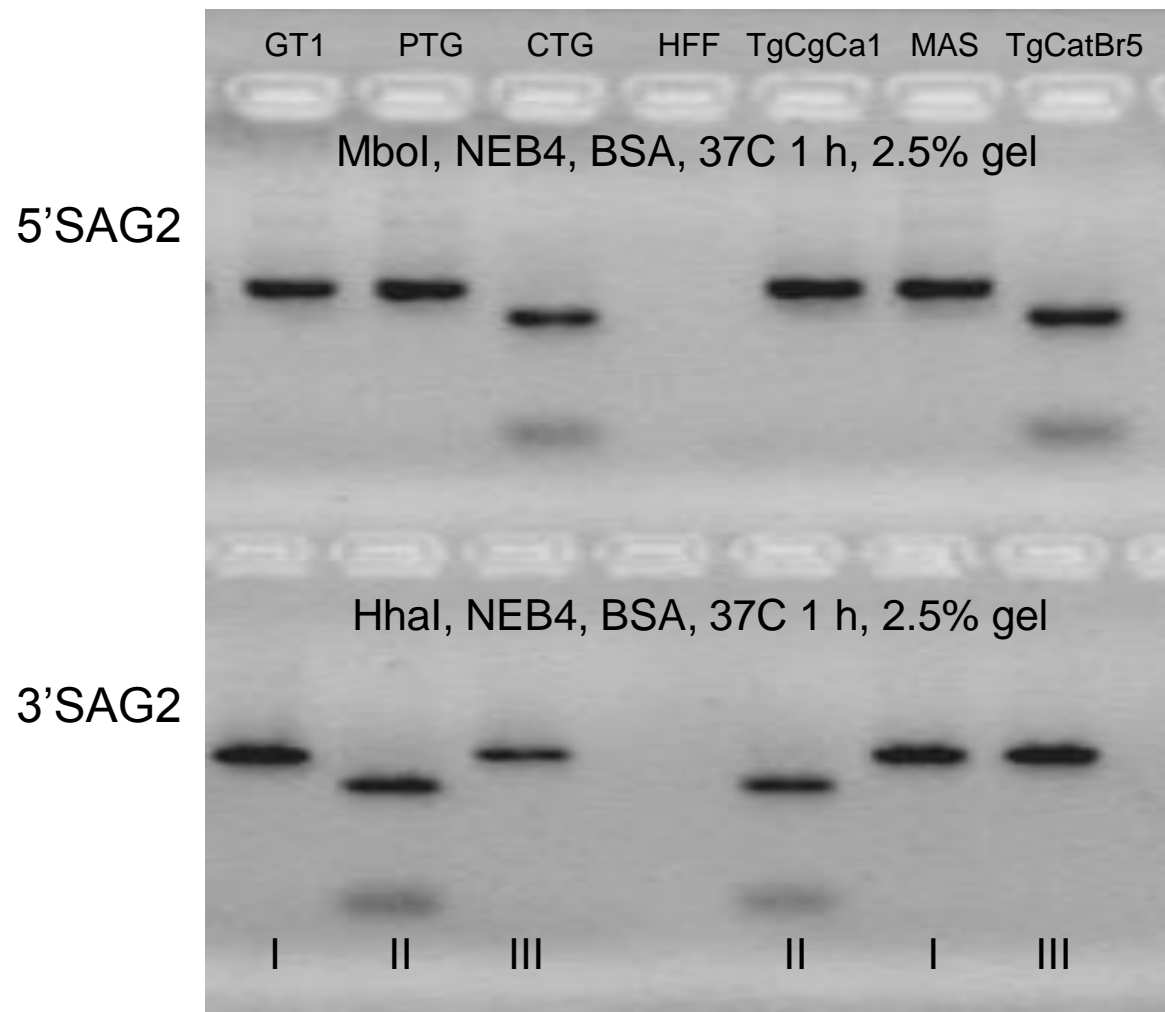


alt. SAG2

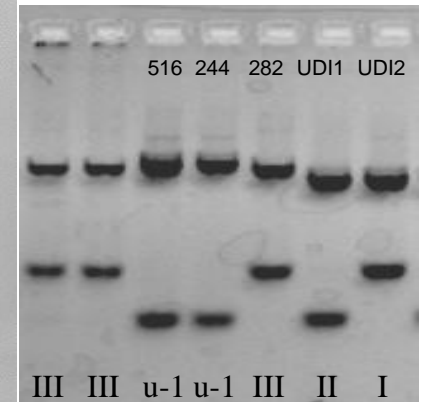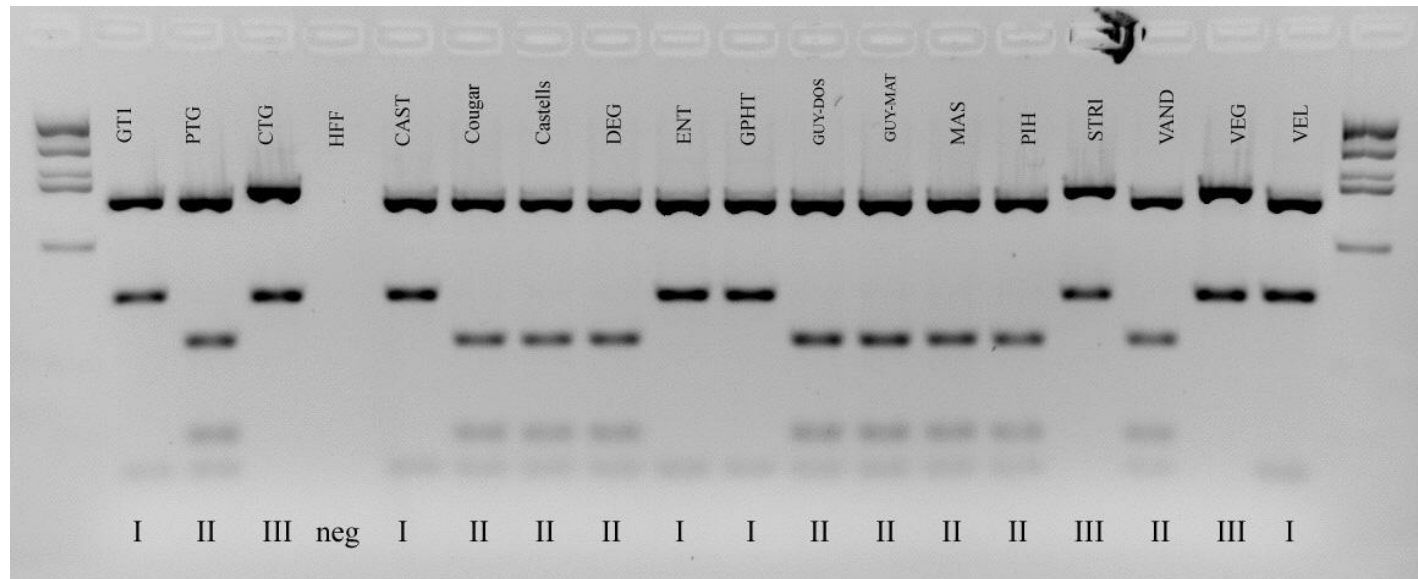

Brazil cats


HinfI + TaqI, NEB3, BSA, 37°C 30 min, 65°C 30 min, 2.5% gel

SAG3

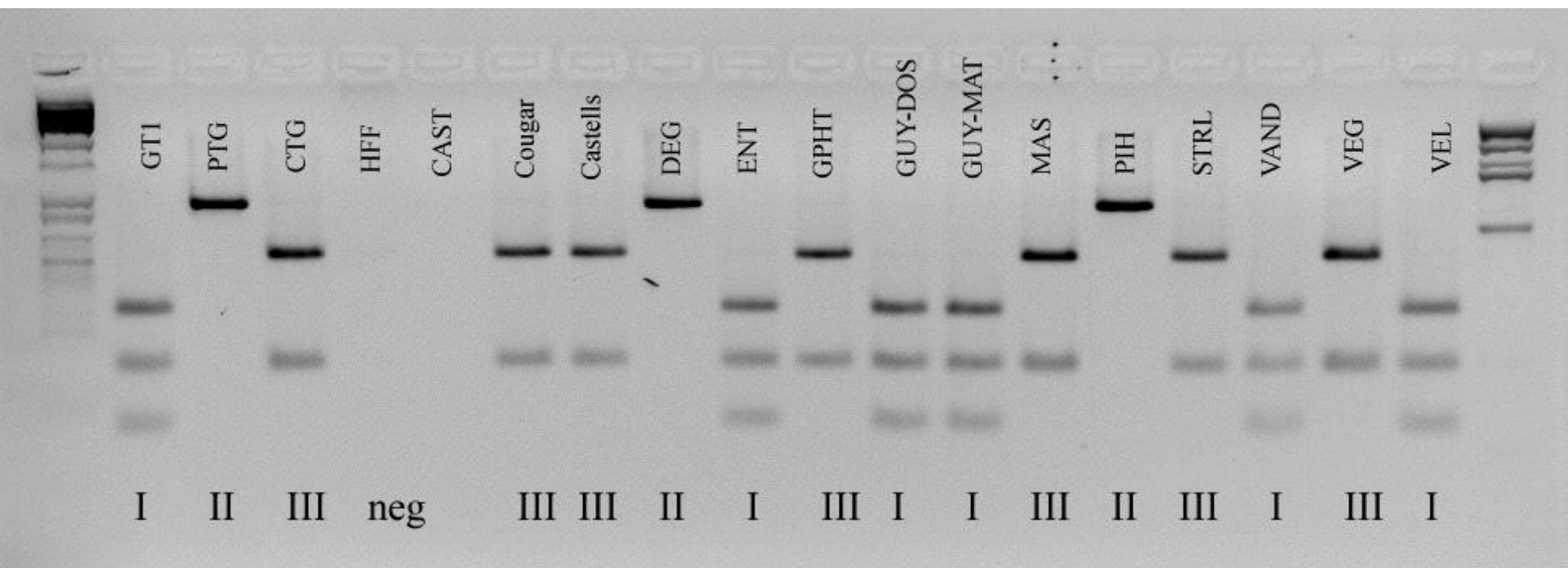


Ncil, NEB4, BSA, 37C 1 h, 2.5% gel

BTUB

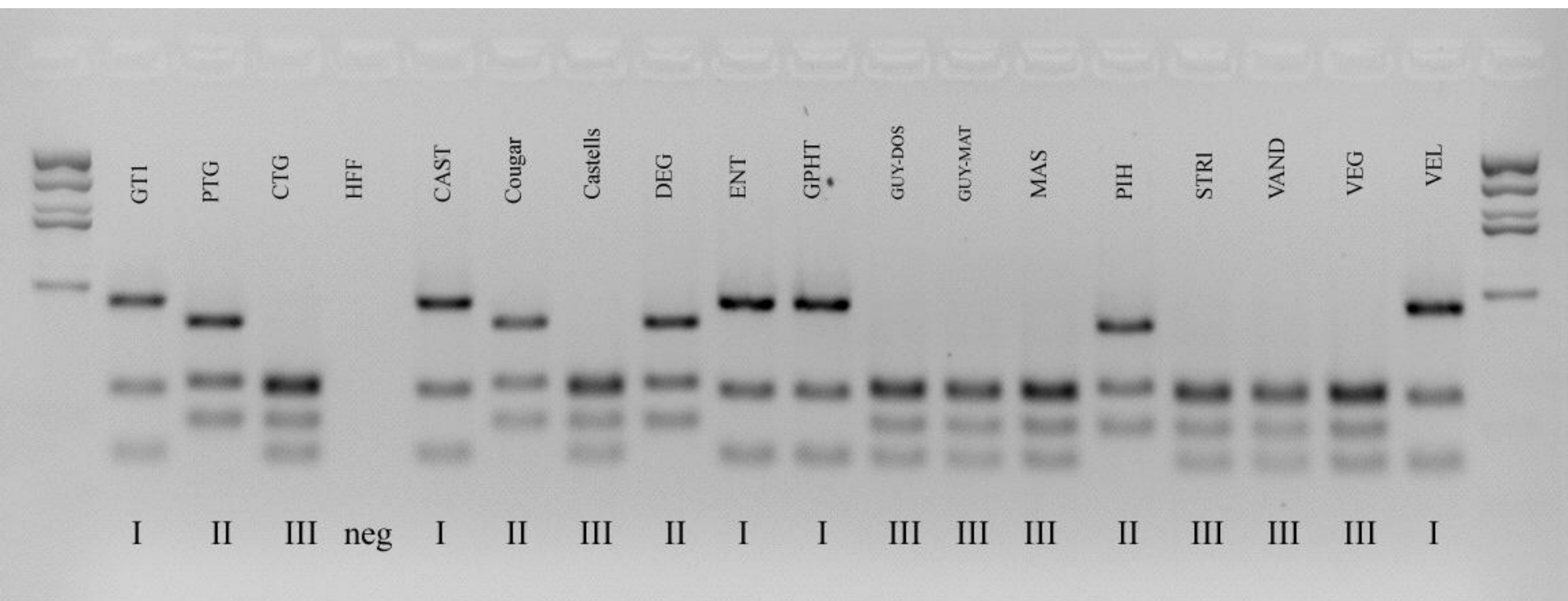

BsiEI + TaqI (double digest), NEB4, BSA, 60C 1 h, 2.5% gel


GRA6

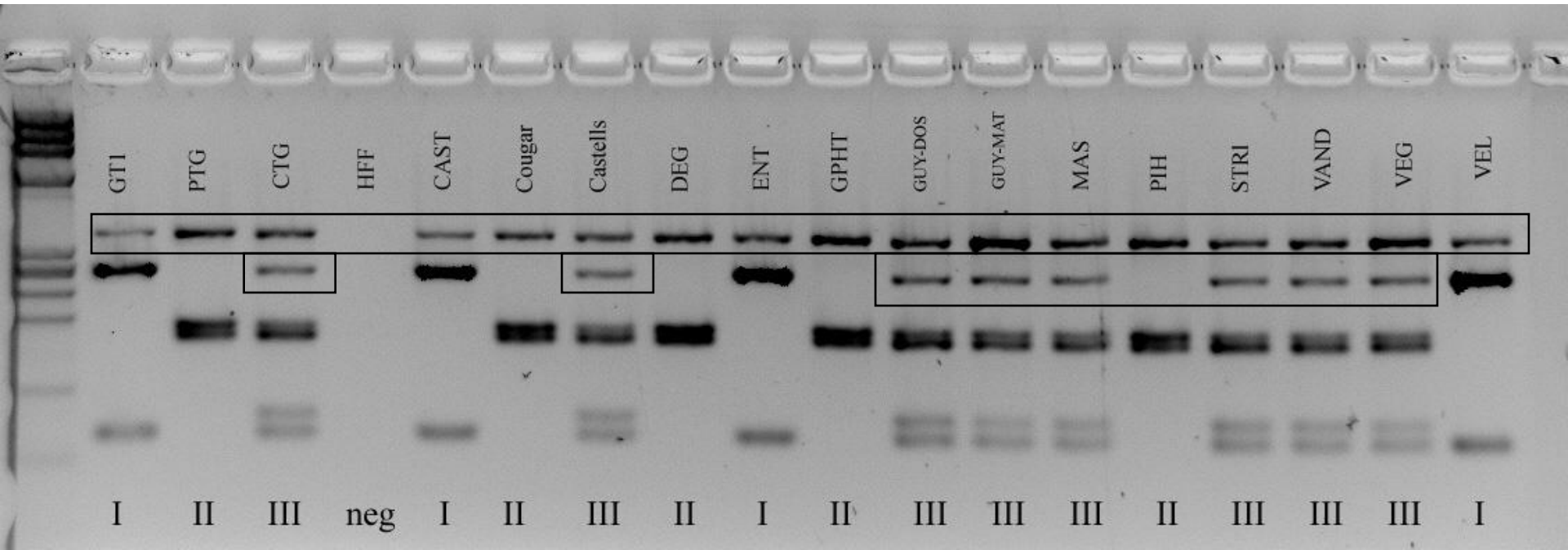

MseI, NEB2, BSA, 37C 1 h, 2.5% gel


Note: Bands in boxes are from incomplete digestion

c22-8

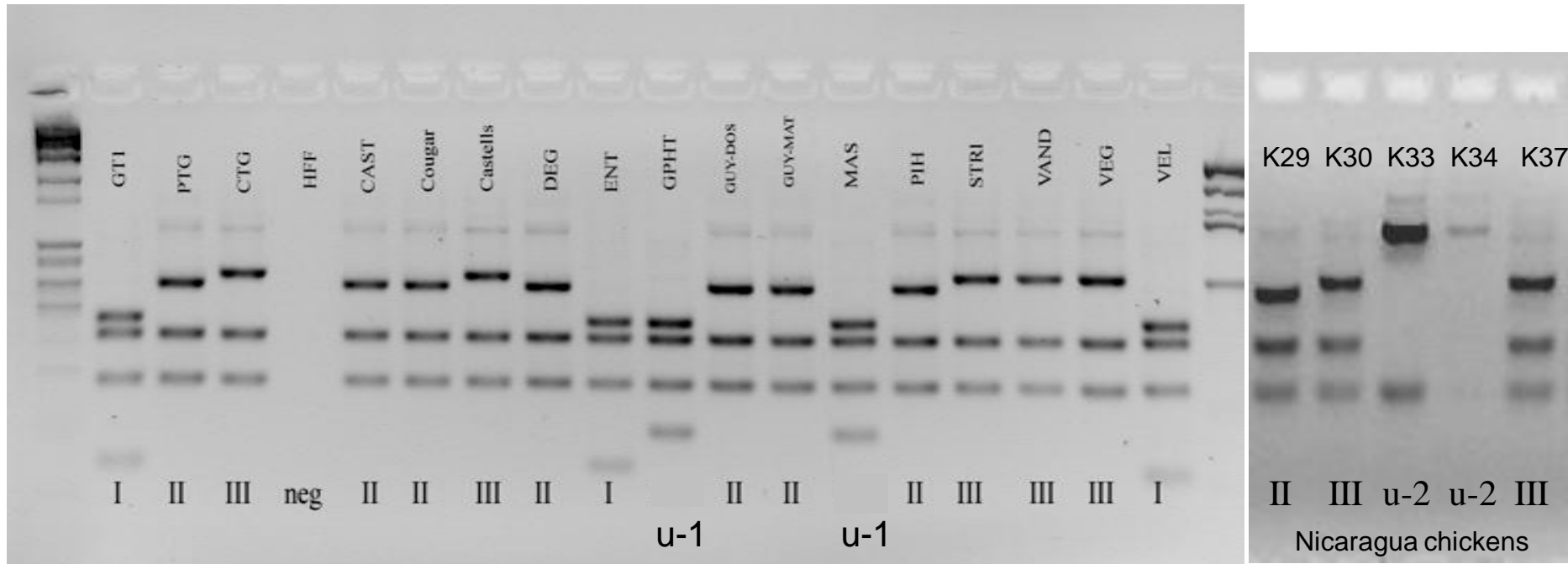


BsmAI + MbolI (double digest), NEB2, BSA, 37°C 30 min, 55°C 30min, 2.5% gel

c29-2

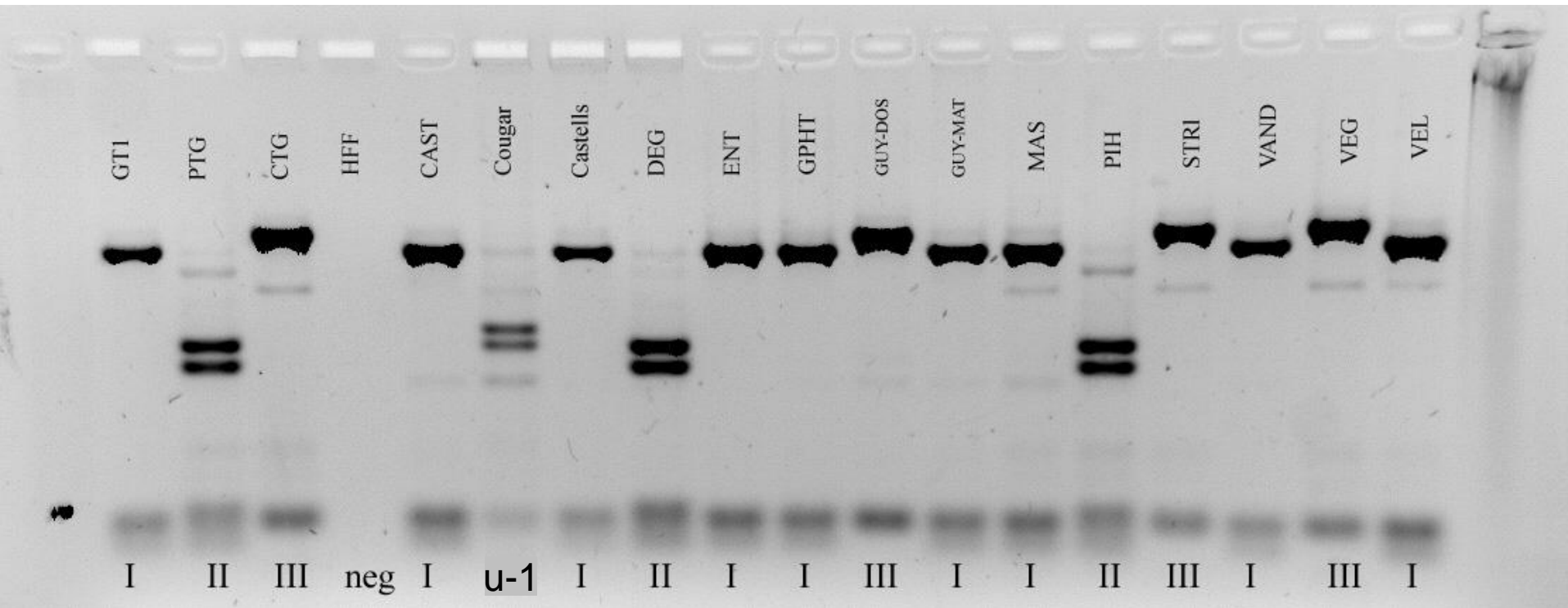


HpyCH4IV + RsaI (double digest), NEB1, BSA, 37°C 1 h, 2% gel

L358

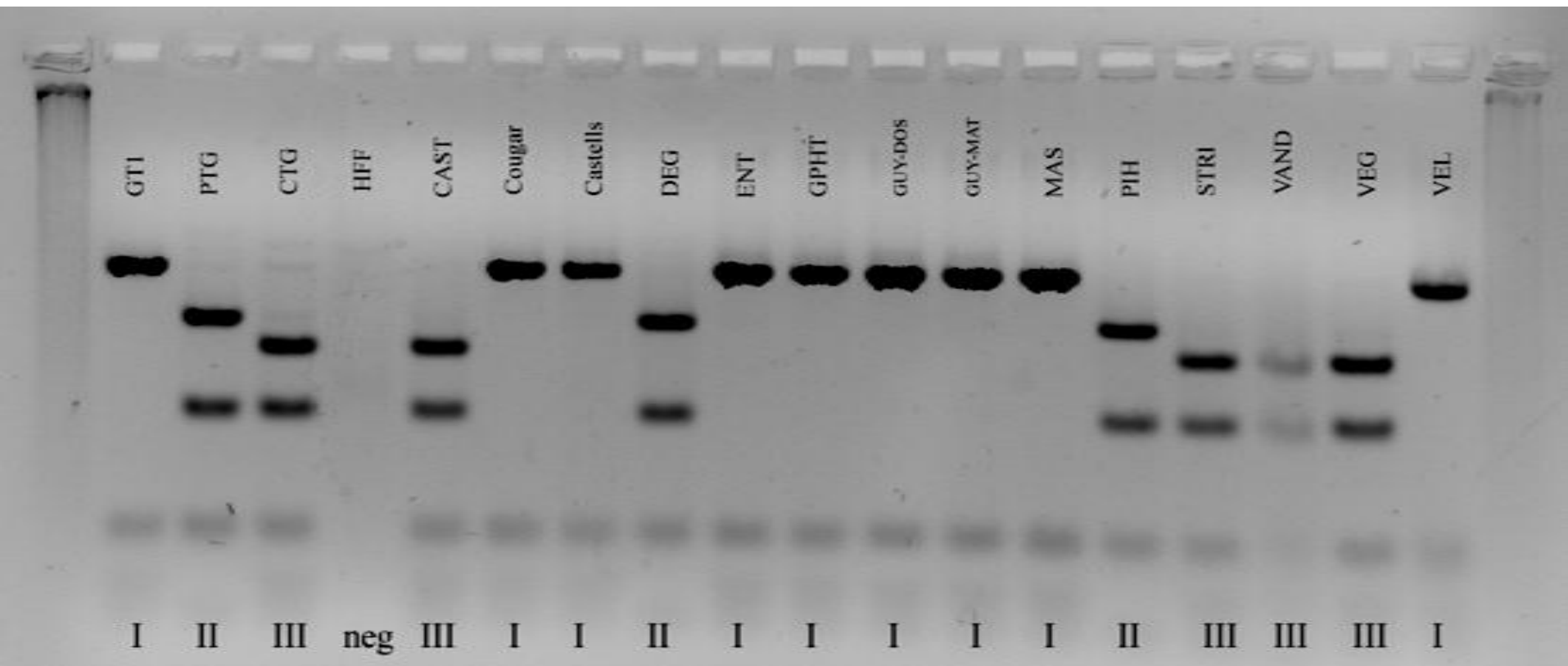


HaeIII + NlaIII (double digest), NEB4, BSA, 37°C 1 h, 2.5% gel


PK1

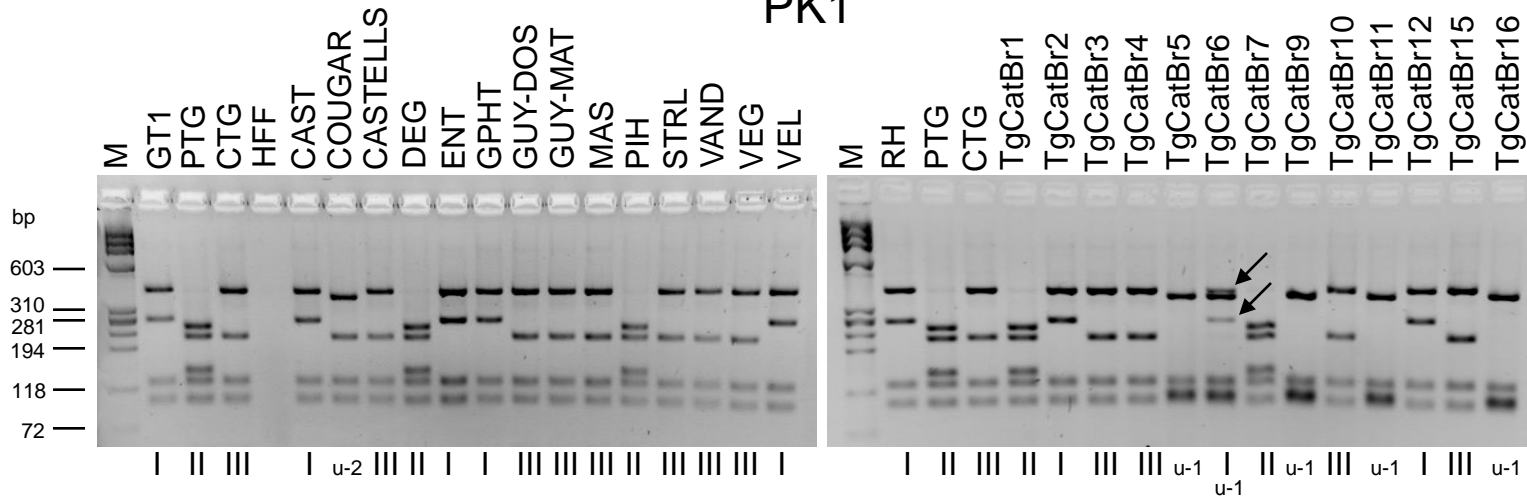

Aval + RsaI (double digest), NEB4, BSA, 37°C 1 h, 2.5% gel

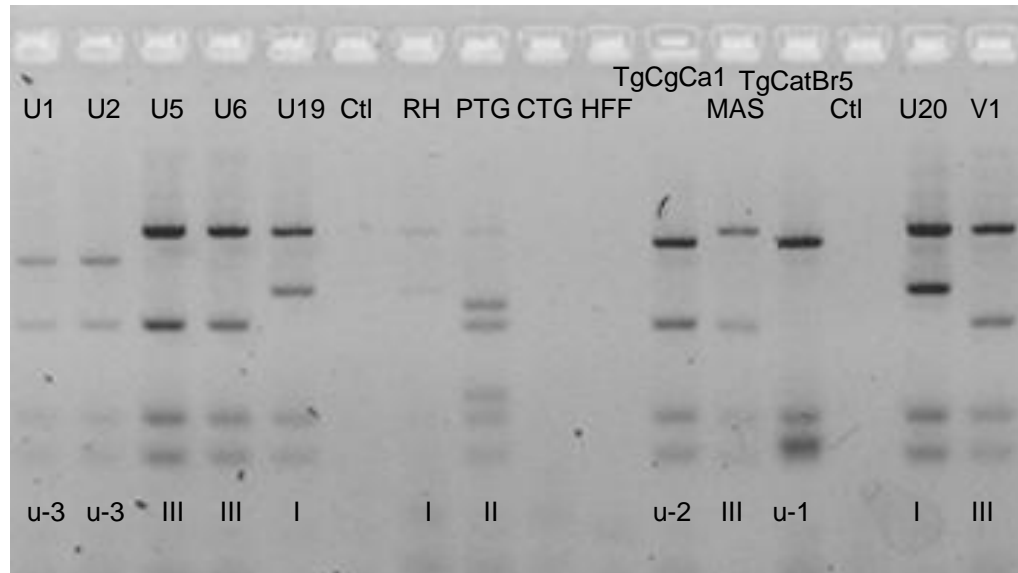

Apico

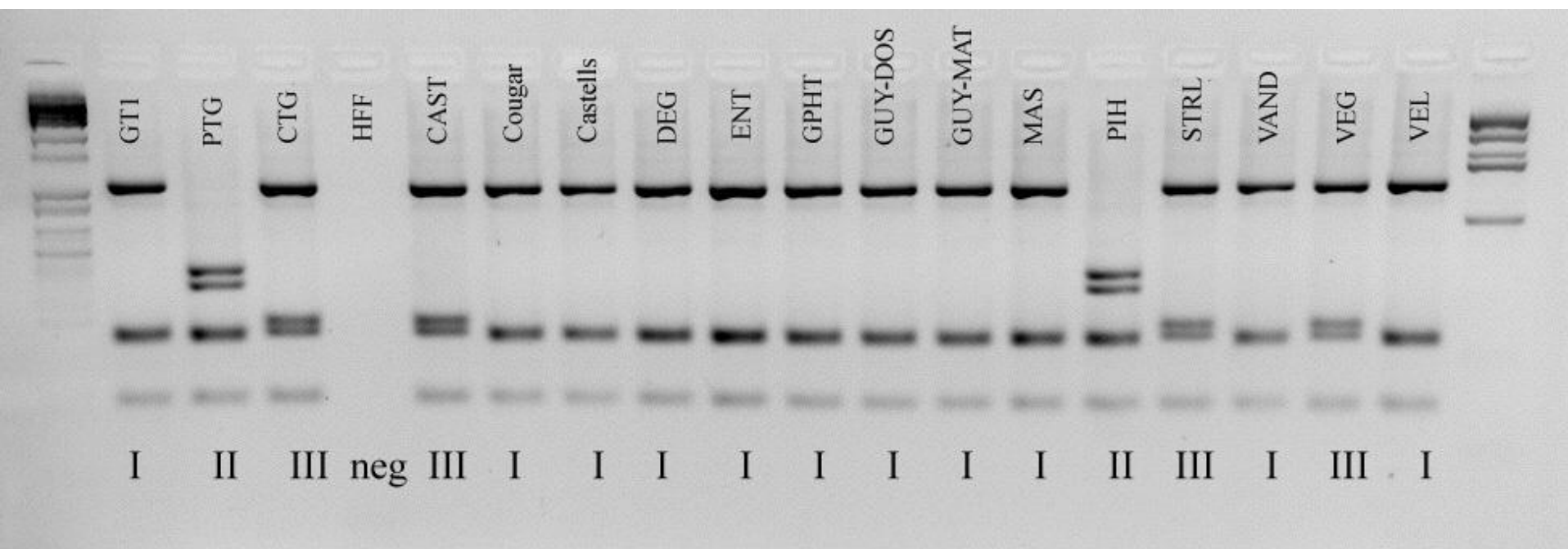

AfIII + DdeI (double digest), NEB2, BSA, 37°C 1 h, 3% gel